



ANTEC DECADE ELITE

Clarity Control Module

ENG

Code/Rev.: M179/90B

Date: 2024-11-13

Phone: +420 251 013 400

clarity@dataapex.com

www.dataapex.com

DataApex Ltd.
Petrzilkova 2583/13
158 00 Prague 5
Czech Republic

Clarity[®], DataApex[®] and ▲[®] are trademarks of DataApex Ltd. Microsoft[®] and Windows[™] are trademarks of Microsoft Corporation.
DataApex reserves the right to make changes to manuals without prior notice. Updated manuals can be downloaded from www.dataapex.com.

Author: MP

Contents

1 Antec DECADE Elite Control Module	1
2 Requirements	2
3 Installation Procedure	3
3.1 Hardware - Wiring	3
3.2 DECADE Elite digital inputs and outputs	3
3.3 Clarity Configuration	6
4 Using the control module	9
4.1 Method Setup - Acquisition	9
4.1.1 Method Setup - Acquisition - Main	10
4.1.2 Method Setup - Acquisition - Output	13
4.1.3 Method Setup - Acquisition - Time Table	14
4.2 Method Setup - Event Table	15
4.3 Method Setup - Advanced	16
4.4 Device Monitor	17
4.5 Antec Detector Setup	19
4.6 Digital Outputs Control	23
5 Report Setup	24
6 Troubleshooting	25
6.1 others.ini File Setup	26

To facilitate the orientation in the **Antec DECADE Elite** manual and **Clarity** chromatography station, different fonts are used throughout the manual. Meanings of these fonts are:

Open File (italics) describes the commands and names of fields in **Clarity**, parameters that can be entered into them or a window or dialog name.

WORK1 (capitals) indicates the name of the file and/or directory.

ACTIVE (capital italics) marks the state of the station or its part.

Chromatogram (blue underlined) marks clickable links referring to related chapters.

The bold text is sometimes also used for important parts of the text and the name of the **Clarity** station. Moreover, some sections are written in format other than normal text. These sections are formatted as follows:

Note: Notifies the reader of relevant information.

Caution: Warns the user of possibly dangerous or very important information.

Marks the problem statement or trouble question.

Description: Presents more detailed information on the problem, describes its causes, etc.

Solution: Marks the response to the question, presents a procedure how to remove it.

1 Antec DECADE Elite Control Module

This manual describes the setting of the **Antec DECADE Elite** and **DECADE Lite** Electrochemical detectors. The control module enables direct control of the instrument over USB or LAN, depending on the communication interface installed in the detector.



Fig. 1: Antec DECADE Elite detector

Direct control means that the detector can be completely controlled from the **Clarity** environment. Instrument method controlling the analysis conditions will be saved in the measured chromatograms.

2 Requirements

- **Clarity** installation with LC Control module (p/n A24).
- Free USB port in the PC or detector connected to the same LAN as the PC.
- LAN cable (p/n SK08) for network-connected detector or USB A-B cable (p/n SK06) for USB-connected one.

Note: Cables are not part of **Clarity** Control Module. It is strongly recommended to order required cables together with the Control Module.

3 Installation Procedure

3.1 Hardware - Wiring

The **Antec DECADE Elite** detector is controlled either by serial (RS232) or LAN communication. For the RS232 communication, there is an embedded USB-to-RS232 port, which requires standard A-B USB cable (p/n SK06) to use. For the LAN communication, standard cross LAN cable (p/n SK08) is used.

3.2 DECADE Elite digital inputs and outputs

The **Antec DECADE Elite** detector is, besides the USB or LAN connector, equipped with several other connectors allowing the use of various digital inputs and digital outputs. The connectors are placed on the back panel of the instrument and shown in the picture below:

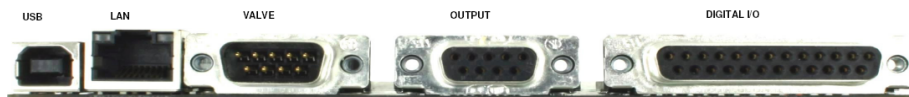


Fig. 2: Back panel of the Antec DECADE Elite detector

There are three input/output connectors installed on the back side of the **Antec DECADE Elite** that are not used for the digital communication. **Valve** is 9-pin DB9M connector, **Output** is a 9-pin DB9F connector and **Digital I/O** is 25-pin DB25F connector. **Valve** connector is used for operating the optionally connected inject valve, **Output** connector hosts analog signal outputs and **Digital I/O** contains pins used for digital inputs and digital outputs of the **Antec DECADE Elite** detector. Both **Valve** and **Output** connectors have the same orientation and organization of pins, but individual pins serve for different purposes.

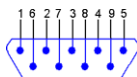


Fig. 3: Valve (DB9M) and Output (DB9F) connector

The meanings of individual pins are as follows:

Tab. 1: Valve connector contacts

Pin No.	Pin Name
1, 5	GND
2	Transmit (TxD)
3	Receive (RxD)
4	-
6	Switch (Handvalve)
7	-

Pin No.	Pin Name
8	-
9	Inject Marker

Tab. 2: Output connector contacts

Pin No.	Pin Name
1	Analog Output 1
2	Analog Output 2
3	Analog Output 3
4	Analog Output 4
5	Analog Output 5
6-9	GND

The Digital I/O connector has different pin allocation:

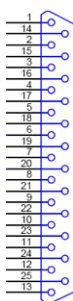


Fig. 4: Digital I/O (DB25F) connector

Tab. 3: Digital I/O connector contacts

Pin No.	Type	Function in Default setting
1	Relay Output	Relay 1 - NC
2	Relay Output	Relay 1 - NO
3	Relay Output	Relay 1 - Common
4	Relay Output	Relay 2 - NC
5	Relay Output	Relay 2 - NO
6	Relay Output	Relay 2 - Common
7	TTL Output	AUX 1
8	TTL Output	AUX 1
9	TTL Output	AUX 1
10	TTL Output	AUX 1
11	TTL Output	Overload
12	TTL Input	Cell 1 On
13	TTL Input	Cell 2 On

Pin No.	Type	Function in Default setting
14	TTL Input	Cell 3 On
15	TTL Input	Cell 1 Off
16	TTL Input	Cell 2 Off
17	TTL Input	Cell 3 Off
18	TTL Input	Cell 1 Autozero
19	TTL Input	Cell 2 Autozero
20	TTL Input	Cell 3 Autozero
21	TTL Input	Cell 1 Start IN
22	TTL Input	Cell 2 Start IN
23	TTL Input	Cell 3 Start IN
24	TTL Input	Cell 4 Start IN
25	GND	GND

Note: Level triggered TTL input contacts require a minimum TTL *LOW* pulse duration of 100 ms. If multiple activations are required, the next pulse should be given after at least 100 ms TTL *HIGH*. When the input is kept *LOW*, only one activation will occur. The default value for the TTL output is *HIGH* (5 V).

The use of digital inputs and outputs in **Clarity** is described in the topic concerning [Antec Detector Setup](#) dialog (for more details see also chapter **Antec Detector Setup** on pg. 19).

Caution: **Clarity** expects that the Inputs and Outputs are set in the default state. Do not change the mapping of inputs and outputs if you want to use them in **Clarity**.

3.3 Clarity Configuration

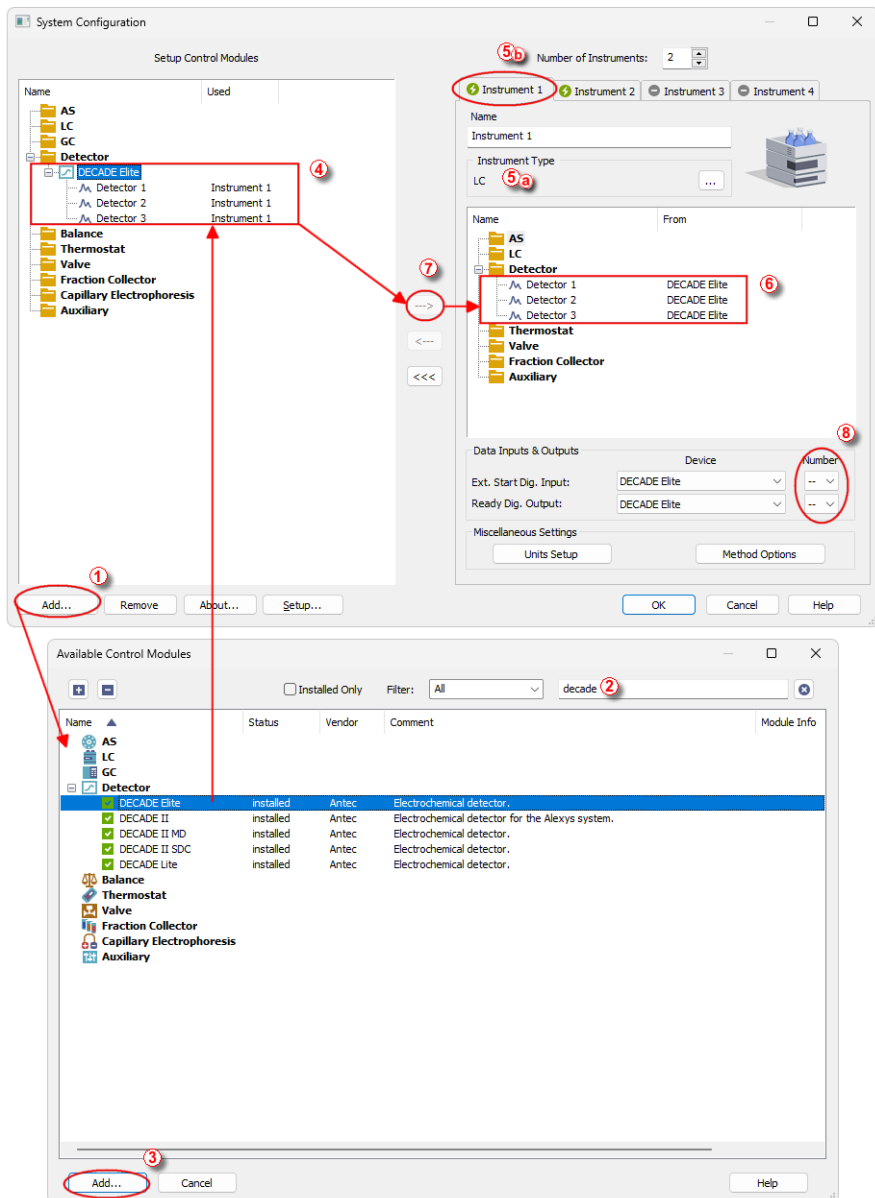



Fig. 5: System Configuration

- Start the **Clarity** station by clicking on the  icon on the desktop.
- Invoke the *System Configuration* dialog accessible from the *Clarity* window using the *System - Configuration...* command.
- Press the *Add* button ① (see 3.3 on pg. 6.) to invoke the *Available Control Modules* dialog.
- You can specify the searching filter ② to simplify the finding of the driver.
- Select the **Antec DECADE Elite** detector and press the *Add* ③ button.

The [Antec Detector Setup](#) dialog will appear.

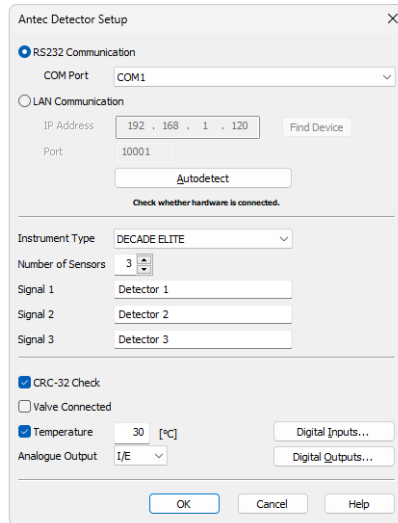
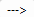


Fig. 6: Antec Detector Setup

- Select the correct type of the *Communication* and fill in the appropriate *COM Port* or *IP Address*, then press the *Autodetect* button. You can also choose your detector from the list after pressing *Find Device* button. If the communication is correct, model type including the number of configured cells will appear below the *Autodetect* button. Press the *OK* button to close the [Antec Detector Setup](#) dialog.

Note: The [Antec Detector Setup](#) dialog is more closely described in the chapter "**Antec Detector Setup**" on pg. 19.

The **Antec DECADE Elite** item will appear in the *Setup Control Modules* list of the *System Configuration* dialog.

- Change the *Instrument Type* ⑤ ③ on the desired *Instrument* tab ⑤ ⑥ to LC and drag the **Antec DECADE Elite** item from the *Setup Control Modules* list on the left ④ to the *Instrument* tab on the right ⑥, or use the  button

⑦ .

- Set the *Ext. Start Dig. Input* and *Ready Dig. Output* numbers ⑧ for your acquisition card according to the wires being used for synchronization.

4 Using the control module

New [Acquisition](#) tabs (one for each detector signal assigned to the Instrument) are created in the *Method Setup* dialog. They can be accessed by switching to the desired signal in the *Select Detector* section on the top of the dialog.

Note: Detector does not support the *Start - Restart* option in the *External Start Options* section which can be found in the lower right part of the *Method Setup - Measurement* tab.

4.1 Method Setup - Acquisition

The *Method Setup - Acquisition* tab for each signal of the **Antec DECADE Elite** detector has three sub-tabs: [Main](#), where basic information on the signal is set, [Output](#) sub-tab with other signal properties and the [Time Table](#) sub-tab containing the *Time table* governing the time program of the detector.

In the right section on all *Method Setup - Acquisition* sub-tabs for **Antec DECADE Elite** detectors two buttons are placed:

From Det

Pressing the button loads all detector parameters as they are set in the **Antec DECADE Elite** detector. These parameters are set in appropriate fields of the *Method Setup - Acquisition* dialog.

Det Status

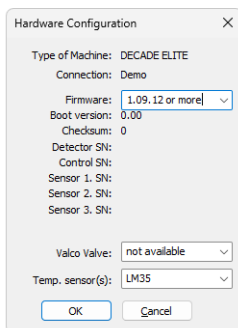


Fig. 7: Hardware Configuration

After pressing the *Det Status* button the *Hardware Configuration* dialog opens. It displays the detector type, communication port through which it is connected, used firmware version, serial numbers of all sensor boards that are installed and other information.

4.1.1 Method Setup - Acquisition - Main

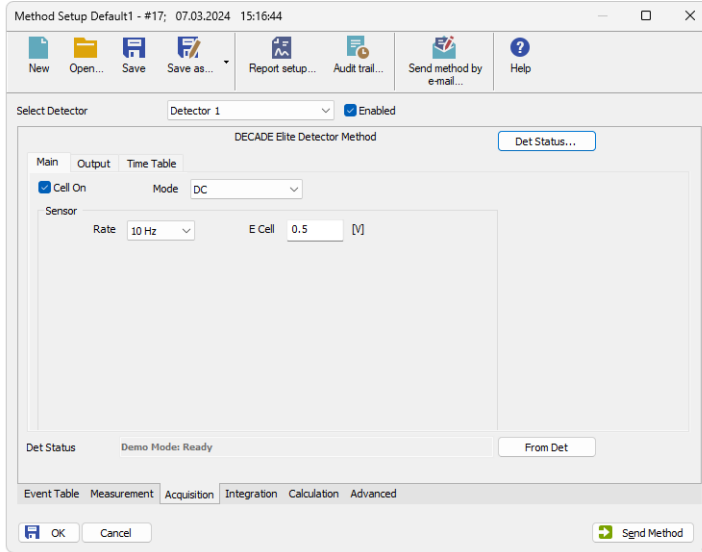


Fig. 8: Method Setup - Acquisition - Main

The *Main* sub-tab of the [Method Setup - Acquisition](#) dialog is used to specify the operation mode of the particular detector and its other parameters.

Cell On

This checkbox enables to switch the particular cell off. The cell is on by default. The analysis cannot be performed while the cell is off, to be able to start it it is necessary to disable the particular detector in the top part of its [Method Setup - Acquisition](#) tab.

Mode

The **Antec DECADE Elite** detector can work in three operation modes - *DC*, *Pulse* and *Pulse 2*. The selection of the mode influences the parameters which can be set in the *Sensor* section.

Sensor

This section differs for the *DC*, *Pulse* and *Pulse 2* modes. For *DC* mode only the sampling *Rate* and *E Cell* potential are defined, for *Pulse* and *Pulse 2* modes the section is used for specifying the full potential profile.

DC Mode

This mode makes the detector to measure on a set potential without changing it.

Rate

Defines the sampling rate of the detector. Possible values differ based on setting of the *Filter* parameter on the [Output](#) tab - in case *Filter* is *OFF* or any filter value other than *Raw* is selected, sample rate is locked to a

proper corresponding sample rate. In case *Raw* option is selected there, possible *Rate* ranges between 1 and 100 Hz.

E Cell

Defines the potential set on the measurement cell.

Pulse mode

This mode allows to set simple rectangular pulse shapes with a single measuring potential and up to 4 different potentials used for renewing electrode surface. The potential shape is as follows:

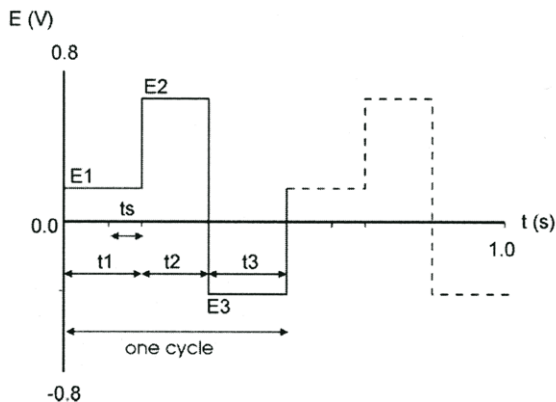


Fig. 9: Pulse mode

Rate

Shows the calculated sample rate of the detector. The calculation is performed as $1/\text{pulse length}$, as calculated from the sum of t_1 , t_2 , t_3 , t_4 and t_5 times.

E1

Defines the measurement potential.

t1

Sets the time for which the measurement potential will be maintained.

E2 .. E5

Sets the potential value for the t_2 .. t_5 time intervals.

t2 .. t5

Defines the duration of the interval for which the E_2 .. E_5 potentials will be set.

Note: The potential set during the t_2 , t_3 , t_4 and t_5 time intervals is used to renew the surface of the electrodes.

ts

Sets the sampling time. The time difference between t_1 and t_s values is used to establish the correct value of E_1 , possible values are always

shorter than $t1$ by at least 60 ms and available in 20 ms increments.

Pulse 2 mode

This mode allows to set more complex pulse shapes with varying measurement times, not necessarily limited to a single potential interval.

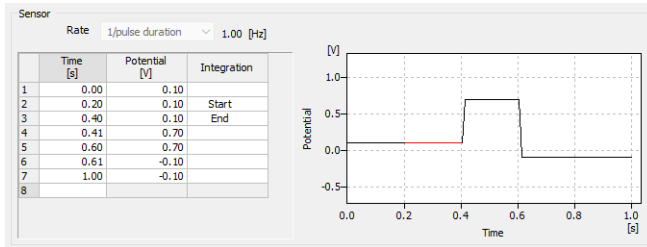


Fig. 10: Pulse 2 mode

Rate

Shows the calculated sample rate of the detector. The calculation is performed as $1/\text{pulse length}$, as received from the last table *Time* value.

Time-Potential table

Allows to set the points where any given potential should be set as pairs of *Time* and *Potential*. The values in the table are interpolated, so in order to obtain as square pulse shapes as possible, use minimum time value differences (0.01 s). Moreover, it is necessary to select the measurement period by selecting its start and end in the *Integration* column. Only a single sampling period per pulse duration is allowed.

Pulse shape graph

Displays the shape of the pulse as set in the time table. The Integration period is displayed in the graph highlighted by red color.

Note: Validity of the pulse shape is dependent on several options. In case the **Antec DECADE Elite** detector drops to the **PAD OVERLOAD** error state after sending the method or start of the measurement, try raising the *Range* parameter on the [Output](#) tab.

4.1.2 Method Setup - Acquisition - Output

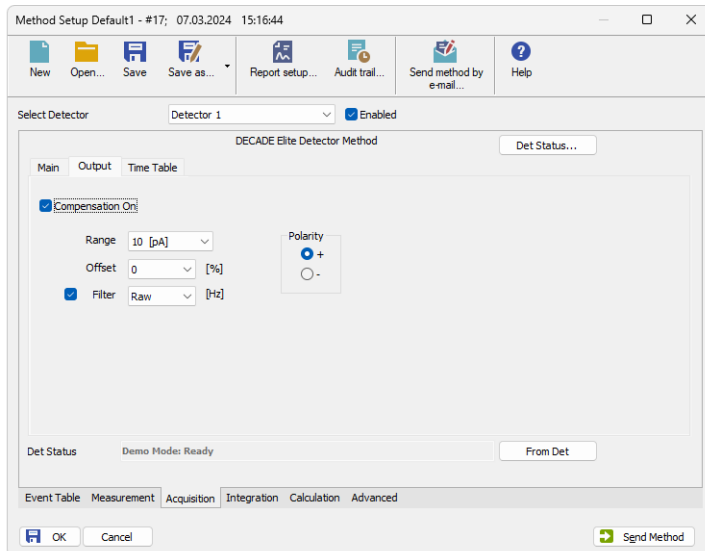


Fig. 11: Method Setup - Acquisition - Output

This sub-tab contains other signal parameters of the particular detector. These are:

Compensation On

This checkbox releases or applies an Autozero offset. When unchecked, any applied autozero offset will be released and the uncompensated cell current will be applied to the output.

Range

Sets the range of the detector. The range can be also changed later as a part of the time program of the detector on the [Method Setup - Acquisition - Time Table](#) sub-tab.

Offset

Sets the y-axis offset of the given signal. This offset is set in percentage of the signal range.

Filter

The *Filter* checkbox allows to use the noise filter for the particular signal. The dropdown serves for specifying the cut-off frequency of the filter (possible values are *RAW* and range between 0.001 and 10 Hz for the *DC* mode and between 0.001 and 0.5 Hz for the *Pulse* mode).

Note: In DC mode the Sensor Rate (Acquisition, Main tab) is fixed for all settings except Filter ON and RAW combination.

Polarity

Sets the polarity of the signal. This parameter may be used for simple switching of negative peaks to positive values.

4.1.3 Method Setup - Acquisition - Time Table

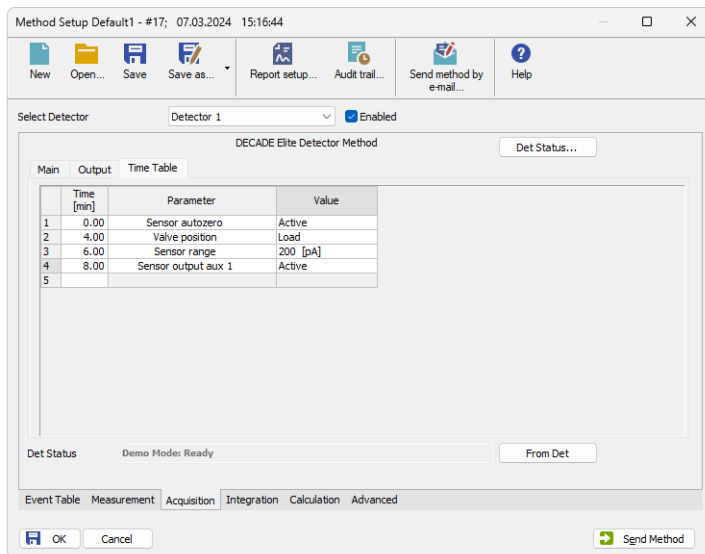


Fig. 12: Method Setup - Acquisition - Time Table

The Time table allows to set timed events on the given detector. The following parameters and outputs can be programmed in the time table:

Valve position, Sensor range, Sensor autozero, Sensor offset, Sensor E cell, Sensor output aux 1(2), Sensor output relay 1(2).

Some events set in the Time table will alter the original values possibly set on the [Method Setup - Acquisition - Main](#) or [Method Setup - Acquisition - Output](#) sub-tabs, others govern the digital outputs of given detector or controller board set in the [Antec Detector Setup](#) dialog.

After entering the desired *Time* for the event (in minutes), it is possible to fill in the *Parameter* field with any possibility from the drop-down list. Add the desired value of the parameter in the *Value* column (if applicable).

4.2 Method Setup - Event Table

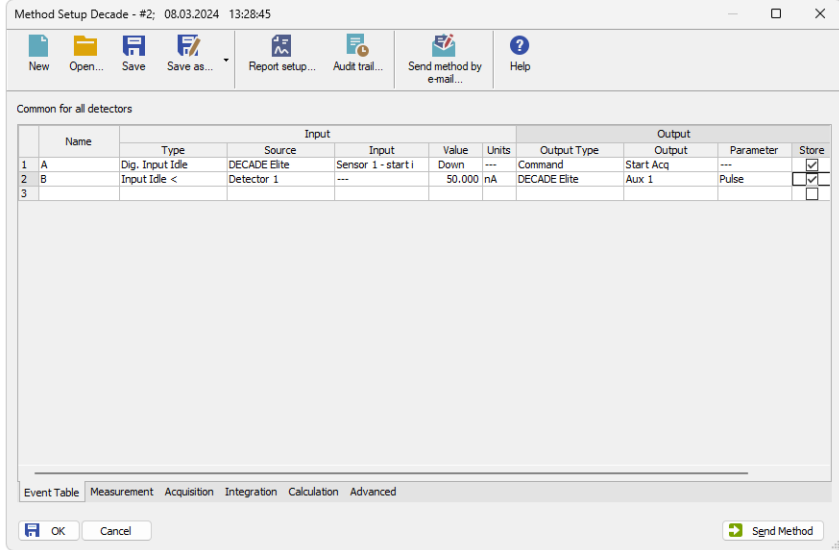


Fig. 13: Method Setup - Event Table

Both the inputs and outputs of the **Antec DECADE Elite** detector can be used in the **Event Table**, either as the input conditions or for the output reaction.

4.3 Method Setup - Advanced

The **Antec DECADE Elite** detector also provides auxiliary signals, which can be added to the resulting chromatograms in the *Method Setup - Advanced* dialog.

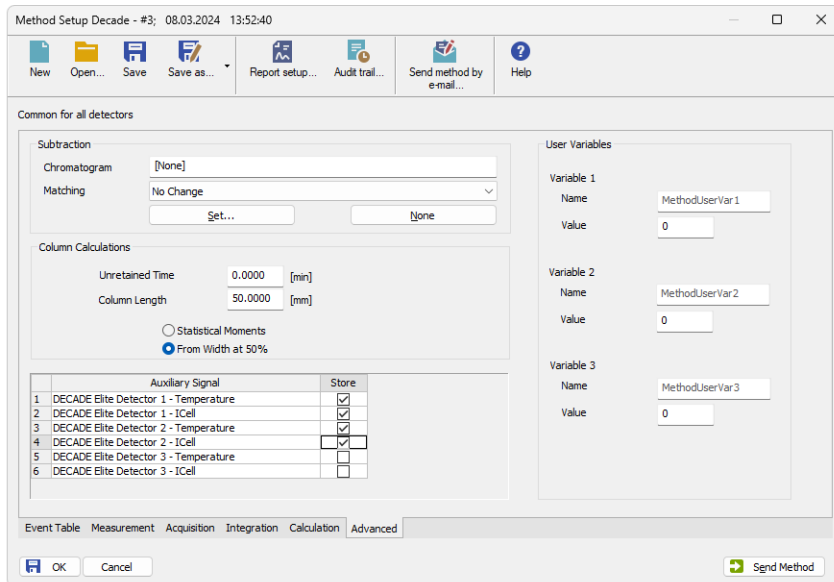



Fig. 14: Method Setup - Advanced

Check the checkbox in the *Store* column to save the particular auxiliary signal to the chromatogram.

4.4 Device Monitor

The window with the detector status can be invoked by the *Monitor - Device Monitor* command from the *Instrument* window or using the  *LC Monitor* icon. It displays the actual signal belonging to all signals configured on the Instrument, as well as other information such as the detector temperatures. Each detector has its own section in the *Device Monitor* window.

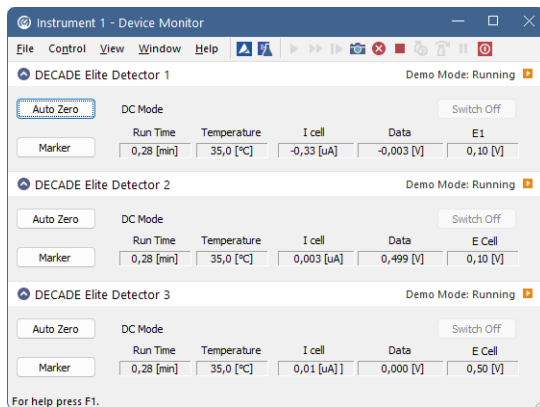


Fig. 15: Device Monitor

The detector can be monitored and partially controlled from the *Device Monitor* window. Available control commands and fields are:

Auto Zero

This button performs the autozero operation - the output signal of the detector is set to zero.

Switch Off/Switch On

This button switches the cell *OFF* or *ON*, depending on the previous status of the cell.

Caution: Clarity requires the cell to provide data, so if the cell is off, button says *Switch On*, it will trigger *NO DATA* error or (in case there is some other detector giving data) *STOPPING* state. Be aware of the fact the External start will start data acquisition even though the cell might be off (button says *Switch On*).

Marker

Invoking this button generates a marker peak in the chromatogram. This function is useful only for analogue plotters as an analysis start marker.

Run Time

Displays the run time of the analysis.

Temperature

Shows the actual temperature in the detector oven of the **Antec DECADE Elite**.

I cell

Shows the cell current of the detector flow cell. The value corresponds to the uncompensated cell current shown on the LCD display of the detector.

Data

Shows the value of the output signal.

E cell

Shows the value of voltage set in the *E cell* field of the [Method Setup - Acquisition - Time Table](#) sub-tab.

This value field is only present in the *DC* mode, it is replaced by the *E1* field in the *Pulse* mode.

E1

Shows the value of voltage set in the *E1* field of the [Method Setup - Acquisition - Time Table](#) sub-tab.

This value field is present in the *Pulse* mode. It is replaced by the *E cell* field in the *DC* mode.

4.5 Antec Detector Setup

Fig. 16: Antec Detector Setup dialog

RS232 Communication

Sets the communication parameters while the detector is communicating over USB.

Note: **Antec DECADE Elite** detector USB port is in fact simulating COM port to be selected in the *Antec Detector Setup* dialog.

COM Port

Selects the communication port.

LAN Communication

Sets the communication parameters while the detector is communicating over LAN.

IP Address

Sets the IP Address to connect to the **Antec DECADE Elite** detector. The detector communicating over LAN will use DHCP server settings, but it is possible to set it to fixed IP Address.

Note: Please check the **Antec DECADE Elite** detector manual for the description of how to set the IP Address of the detector.

Port

Sets the communication port for LAN communication. The default communication port is *10001*.

Find Device

Opens a dialog for choosing the IP address from the list of all **Antec DECADE Elite** detector connected to the same network.

Autodetect

This button is used to test the connection to the **Antec DECADE Elite** detector using the previously set communication parameters (either RS232 or LAN). Obtained data gained from the auto-detection are then set in the dialog.

Instrument Type

Shows the type of the connected **Antec DECADE Elite** detector.

Number of Sensors

Displays the number of sensor boards (detectors) that are present in the **Antec DECADE Elite** detector.

Signal 1 (.5)

Fields belonging to particular detectors enable to set the relevant detector names.

CRC-32 Check

Enables the CRC-32 checking of incoming data.

Valve Connected

Enables to control the connected valve.

Temperature

Sets the operation temperature of the **Antec DECADE Elite** detector oven.

Analogue Output

Sets the mode how the **Antec DECADE Elite** detector prepares analog output data. The analog output with the scale of 1V may either be created from the raw data before they are processed by the inbuilt A/D converter (option *I/E*) or after they are processed by the A/D converter (option *DAC*).

Digital Inputs...

Displays the *Digital Inputs* dialog which sets the events performed after particular digital inputs of the **Antec DECADE Elite** detector are invoked. The **Antec DECADE Elite** supports up to 13 inputs (depending on the number of individual sensor boards in the instrument), of which any 8 can be used in the **Clarity** station.

Caution: It is highly recommended to set all unnecessary inputs and outputs to the *Not available* state, as it will significantly speed up the communication.

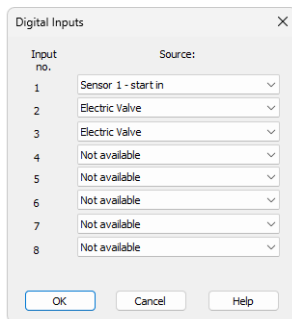


Fig. 17: Digital Inputs

Possible events are:

Not available - digital input with this number is not used.

Manual Valve - option for the setting up of a manual injection valve. The manual injection valve has to be connected to the [Valve Connector](#) located on the back side of **Antec DECADE Elite**. When correctly set in **Clarity System Configuration** the manual injection valve can be used for triggering **Clarity** run if *Ext. Start Dig. Input* set with the same number as Input No. for the manual valve, see ⑧ in **Fig. 5** on pg. 6.. Note that connected manual injection valve can be utilized also through [Event Table](#) after selecting of an option for Digital Input from Type column.

Electric Valve - option for the setting up of electric injection valve. The electric injection valve has to be connected to the [Valve Connector](#) located on the back side of **Antec DECADE Elite** via appropriate RS232 cable. Note the electric valve has to be connected to **Antec DECADE Elite** and switched on before switching **Antec DECADE Elite** on. There are supported two position electric actuators only. When correctly set in **Clarity System Configuration** the electric injection valve can be used for triggering **Clarity** run if *Ext. Start Dig. Input* set with the same number as Input No. for the electric valve, see ⑧ in **Fig. 5** on pg. 6.. Note that connected electric injection valve can be controlled also through [Event Table](#) after selecting of an option for Digital Input from Type column.

Sensor 1(..3) - cell on - monitors the signal which turns the corresponding flow cell on.

Sensor 1(..3) - cell off - monitors the signal which turns the corresponding flow cell off.

Sensor 1(..3) - autozero - monitors the TTL signal which sets output signal of the particular cell to 0.

Sensor 1(..4) - start in - monitors the signal which starts a new time file of the particular detector. **Clarity** doesn't use time files.

Note: To be able to start the analysis in **Clarity** when this signal occurs, it is necessary to set the same signal as *Start IN* signal in the **System Configuration** dialog.

Digital Outputs...

Displays the *Digital Outputs* dialog which sets the events triggered by particular signals of the **Antec DECADE Elite** detector. The **Antec DECADE Elite** supports up to 6 outputs.

Caution: It is highly recommended to set all unnecessary inputs and outputs to the *Not available* state, as it will significantly speed up the communication.

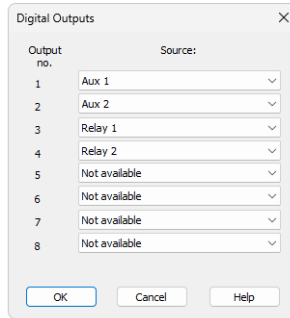


Fig. 18: Digital Outputs

Possible events are:

Not available - digital output with this number is not used.

Sensor 1(..5) - aux 1(2) - output providing a free TTL output contact. Two such outputs are available on each detector board.

Sensor 1(..5) - relay 1(2) - output functioning as a relay contact. Two such outputs are available on each detector board.

4.6 Digital Outputs Control

The **Antec DECADE Elite** detector contains six digital TTL outputs from which the first two are also designed as relay contacts. To set the desired events to actual outputs on the back of the **Antec DECADE Elite** detector use the [Digital Outputs](#) dialog described on pg. 22.

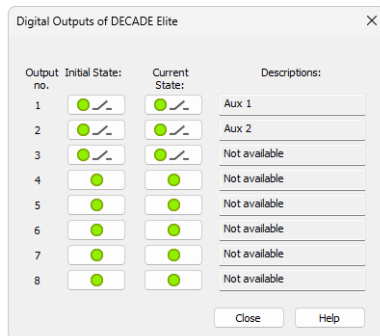


Fig. 19: Digital Output Control

This dialog allows to change the *Initial* and *Actual* state of the digital outputs of the **Antec DECADE Elite**. The digital outputs may be also controlled from the *Method Setup - Event Table* tab.

5 Report Setup

All of the detector settings accessible on the [Method Setup - Acquisition](#) tab for the given signal are reported. To do so, the *Instrument Control* parameter on the *Method* tab of the *Report Setup* dialog must be checked.

The screenshot shows a 'Print Preview' window with a toolbar containing 'Print', 'Print to PDF', 'Send PDF', and 'Close'. The main content area displays detector settings:

Note :
 Autostop : None External Start : Start Only, Down
 Detector 1 : Detector 1 Range 1 : Bipolar, 100 nA, 10 Samp. per Sec.
 Detector 2 : Detector 2 Range 2 : Bipolar, 10000 nA, 10 Samp. per Sec.
 Detector 3 : Detector 3 Range 3 : Bipolar, 10000 nA, 10 Samp. per Sec.
 Detector 4 : Detector 4 Range 4 : Bipolar, 10000 nA, 10 Samp. per Sec.
 Detector 5 : Detector 5 Range 5 : Bipolar, 10000 nA, 10 Samp. per Sec.
 Subtraction Chromatogram : (None) Matching : No Change

Time and Input Events Table (Antec: DECADE)

Name	Type	Input				Output	
		Source	Input	Value	Units	Output Type	Output
A	Dig. Input Idle	DECADE Elite	Sensor 1 - start in	Down	---	Command	Start Acq
B	Inout Idle <	Detector 1	---	50,000	nA	DECADE Elite	Aux 1

Page 1

Fig. 20: Report Setup

6 Troubleshooting

When the remedy for some problem cannot be discovered easily, the recording of communication between **Clarity** and the detector can significantly help the **DataApex** support to discover the cause of the problem.

The recording can be enabled by adding or amending the COMMDRV.INI file in the **Clarity** installation directory (C:\CLARITY\CFG by default). The file can be edited in any text editor (e.g. Notepad). One of the following section, based on the type of communication, should be edited or added:

LAN communication:

```
[TCP_IP 192.168.111.63:10001]
echo=ON
textmode=ON
filename=DecadeElite_%D.txt
reset=OFF
```

RS232 communication:

```
[COM1]
echo=ON
textmode=ON
filename=DecadeElite_%D.txt
reset=OFF
```

Note: On the first line, you may need to correct the IP address or COM port used to communicate with the **Antec DECADE Elite** detector. This number can be displayed in the [Method Setup - Acquisition](#) dialog with *Det Status...* button.

Note: %D (or %d) in the filename parameter means that the log will be created separately for each day. The *reset=off* parameter disables deleting the content of the log each time the station is started during the same day.

The created *.TXT files will greatly help in diagnosis of unrecognized errors and problems in communication.

6.1 OTHERS.INI File Setup

In the OTHERS.INI file (by default located in the C:\CLARITY\CFG directory) user can create section [DecadeII] containing keys, which are influencing the way Clarity deals with missing timestamps in incoming data. The section can contain following keys:

timestamps

Defines what action the control module will do in case missing timestamps problem appears.

Values:

0 - to ignore missing timestamps (Default)

1 - to record missing timestamps into chromatogram audit trail

2 - to calculate and record missing timestamps

start

Defines when the module will start do action from "timestamps"

Values:

0 - Starts from first timestamp, and first expected timestamp is 1 (Default)

1 - the first timestamp will store and starts from second timestamp

2..X - the X timestamp will store and starts from X+1 timestamp

length

Defines maximal count of missing timestamps, which are considered as error

Values:

0 – any number of missing timestamps are considered as error (Default)

1..X - number of missing timestamps exceeding this value will not be considered as error

Example of section (can be directly copied/pasted to OTHERS.INI):

[DecadeII]

; What to do with missing timestamps

**; 0 - Ignore, 1 - Record in audit trail, 2 - Calculate and record
timestamps = 0**

; When to start executing action set in previous key

; 0 - From first missing, X - From X+1 missing

start = 0

; Maximal count of missing timestamps considered as error
; 0 - Any number, X - number of missing timestamps exceeding this value
will not be considered as error
length = 0